*Jade*  
with heating

< relax in style

**QUATTROMED IV-S**

## < relax in style

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## < relax in style

Since being established in 2000, the name CASADA stands for high quality products on the health-, spa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide.

Constant further development, an unmistakable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

## *CONGRATULATIONS!*

---

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Quattromed IV-S.

## *TABLE OF CONTENTS*

---

12	Safety instructions
12	Start up
13	Technical data
14	Equipment / Functions
15	Remote control
16	Massage techniques
16	Troubleshooting
17	Care instructions
17	Storage
17	Conformity declaration
17	Warranty conditions

## SAFETY INSTRUCTIONS

Please carefully read the following instructions to before using this fitness equipment to ensure its proper function and optimal effect. Please retain these instructions for use!

- This massage device complies with the recognised principles of technology and current safety regulations.
- This equipment is maintenance free. Any necessary repairs must be performed by an authorised professional. Improper use and unauthorised repairs are prohibited for safety reasons and will void the warranty.
- Never touch the power plug with wet hands.
- Please avoid the equipment coming into contact with water, high temperatures and direct sunlight.
- Be sure to unplug the equipment from the power point during extended periods of non-use to eliminate the risk of a short circuit.
- Do not use damaged plugs, cords or loose sockets. Damaged plug or cords must be replaced by the manufacturer, service representative or qualified personnel. In the event of malfunction unplug the equipment immediately. We assume no liability for damages resulting from wrongful or improper use.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: head diseases recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.
- Do not use the device while driving.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.

## START UP

Please connect the adapter with the device and subsequently plug the power plug into the power socket. Now turn on the device by the remote control.

## TECHNICAL DATA

**Measurements:** 115 x 44 x 12 cm

**Weight:** 9 kg

**Power source:** Input: AC 100-240 V ~ 50/60 Hz  
Output: DC 24 V, 2000 mA

**Rated power:** 25 watt

**Operating time automatic:** 15 minutes

**Certificates:**

## EQUIPMENT / FUNCTIONS

Please position the Quattromed IV-S on a stable chair. The burden of your bodyweight intensifies the massage effect. Please ensure that your whole body is positioned centrally so that the massage heads run parallel to your spine.

The Quattromed IV-S has an S-Line mechanic, which is based on the natural shape of the spine. This especially benefits your lumbar region by massage heads, which are particularly close to your body. The natural forces of jade stones supports your positive effect of this massage device. By the integrated heat function whose effect is intensified by natural stones, you feel a particularly soothing, warm and comfortable massage.


The Quattromed IV-S has a detailed amorous design that looks great anywhere. Say goodbye to hardened muscles and let yourself be spoiled by your Quattromed IV-S.

The massage function will switch off automatically after 15 minutes. In order to increase the lifespan of the motors, please switch on the massage heads again only after a minimum of 30 minutes have passed. This gives the motors time to cool down again. You should also allow a period of rest for your body. We recommend conducting the massages only two, maximum three times a day in order to not overstress your muscles.

For your own safety we equip the machines with a protection against overheating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk, the machines will switch off automatically and will only allow switching on again after an appropriate downtime between 30 and 40 minutes.

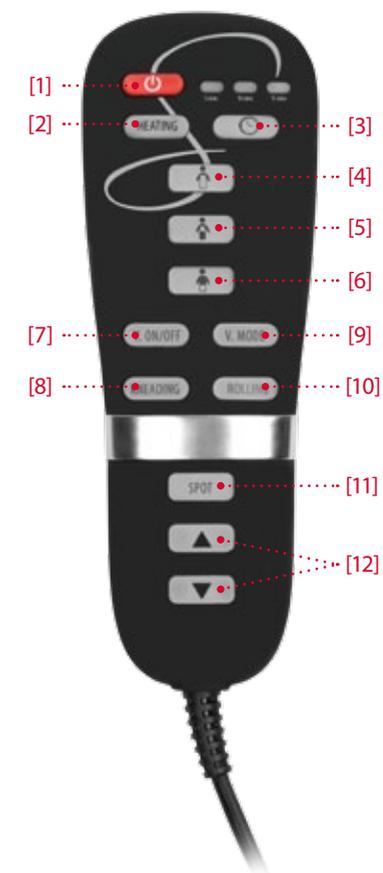**[1] On/Off**  
Switch the device on/off here. When switched off, the massage stops and the mechanism returns automatically to the lower rest position.

**[2] Heating**  
Switch the heating function on/off.

**[3] Time (5-15 minutes)**  
The preset massage time is 15 minutes. This button sets the desired massage duration. (5 minute intervals / 5-10-15).

#### SELECT MASSAGE AREA

**[4]** Switch the massage on/off over the entire back area.

**[5]** Switch the massage on/off in the shoulder and neck area.

**[6]** Switch the massage on/off in the back and lumbar area.

**[7] V. On/Off**  
Switch the vibration massage in the seat area on/off.

**[8] Kneading**  
This button switches the kneading function on/off.

**[9] V. Mode**  
This button controls the intensity of the vibration. (gentle - medium - strong)

**[10] Rolling**  
This button switches the rolling function on/off.

**[11] Spot**  
Press this button once to activate a spot massage.

**[12] ▲/▼**  
These buttons position the massage rollers in the desired section of the back during a spot massage. By holding down the buttons the massage rollers move upwards or downwards. Once the rollers have reached the desired position, let go of the button.

## MASSAGE TECHNIQUES

### SHIATSU/KNEADING

Shiatsu (finger pressure) is a form of body therapy developed in Japan. The massage technique consists of gentle rhythmic twists and rotations with deep effect. Shiatsu aims to stimulate individual points in the body and to mobilise the muscles.

### VIBRATION

The Vibration puts entire body regions into even vibrations, which continue all the way into the depths. Stimulating for circulation, warming and relaxing for the muscles.

### ROLLING

The ROLLING is a pleasant change between tension and release of the muscles and therefore is especially calming and relaxing. This technique is often used as a conclusion of the massage for relaxation of the muscles and for regeneration.

### HEATING

A deep and warmth-radiating effect, where one can relax perfectly. This function is especially enjoyable in case of muscle tension due to its deep effect.

## TROUBLESHOOTING

### Error

The massage device cannot operate even though it has been connected to the power supply.

### Adjustment

- Check the connection of the adapter.
- The device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and allow the device to cool off for 30 to 50 minutes.

## STORAGE / CARE

Regular care will extend the life of your Quattromed IV-S whilst also ensuring safe operation at all times.

1. Recommended storage: dry, free of dust.
2. Keep liquids and open fire away from the Quattromed IV-S
3. Clean the equipment using a dry, soft cloth.
4. Always unplug the Quattromed IV-S when not in use.

## CONFORMITY DECLARATION

Corresponds to the European guidelines 2006/95/EC.

## WARRANTY CONDITIONS

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer. The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Casada International GmbH  
Obermeiers Feld 3  
33104 Paderborn, Germany  
phone (+49) (0 52 54) 9 30 98 - 0  
fax (+49) (0 52 54) 9 30 98 - 30  
info@casada.com  
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright (c) 2014 Casada International GmbH. All rights reserved.

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